

ReSound Live™

Product Information

Product Description

ReSound Live gives your hearing a surround sound experience. It's like going from a set of stereo speakers to a full blown surround sound system.

The surround sound experience is a greater sense of awareness – of 'being there'. Surround sound provides a more full and detailed sound quality, increased ability to locate where sounds are coming from and the ability to hear every word clearly – even in challenging environments.

Careful matching of ReSound Live to mimic the natural human ear is a key element in recreating the surround sound experience.

ReSound Live has two state-of-the-art directionality systems that work together with the surround sound processor. The new AutoScope Adaptive Directionality™ and Natural Directionality™ II set new standards for conventional and asymmetrical fittings.

The powerful Dual Stabilizer® II DFS with WhistleControl™ ensures unprecedented distortion free feedback suppression even in dynamic environments. Combined with the most robust noise reduction system (NoiseTracker™ II and effective suppression of wind noise) the wearer will get the ultimate in sound quality.

Sound by Design: Small but powerful, the new CIC with remote mic offers both increased invisibility, comfort & power.



LV10-M / LV710-M / LV510-M
LV10-MP / LV710-MP / LV510-MP
CIC


	Live 9	Live 7	Live 5
Personalised blending point - adjustable*	●		
Surround sound processor with integrated wind noise suppression	●	●	●
NoiseTracker™ II noise reduction levels (-3, -6, -8, -10 dB)	●	⊙	○
17-band Warp™ compression with extended bandwidth	●	●	
9-band Warp™ compression with extended bandwidth			⊙
WhistleControl™ – sensitivity levels	●	⊙	○
Dual Stabilizer® II DFS	●	●	●
Impulse Noise Smoother	●	●	●
Natural Directionality™ II*	●		
AutoScope Adaptive Directionality™*	●		
MultiScope Adaptive Directionality™*	●	⊙	
Adaptive Directionality*	●	●	●
Fixed Directionality*	●	●	●
SoftSwitching™*	●	●	●
EchoStop™*	●	●	●
Environmental Optimizer™	●		
Onboard Analyzer™ II data logging	●	●	●
Gain handles	9	7	6
Fully flexible programs**	●	⊙	○
SmartStart™	●	●	●
Acoustic indications of user control and notifications	●	●	●
Complete family of appealing, discreet and comfortable design	●	●	●

○ Standard
⊙ Advanced
● Ultimate


* Not available in single mic devices

** Excluding the CIC devices that do not contain a push button

Fitting Range LV10-M

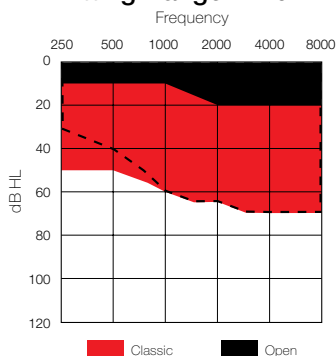

Fitting Range LV10-MP

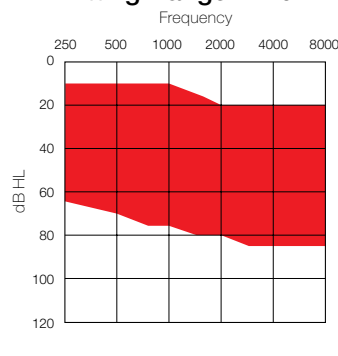

GN ReSound A/S
Lautrupbjerg 7
DK-2750 Ballerup, Denmark
Tel.: +45 45 75 11 11
Fax: +45 45 75 11 19
www.resound.com

Standard Configuration

- Size 10A battery
- Battery door with integrated on/off switch
- Supports StepVent
- Available in 3 colours

Fitting Requirements

- Aventa™ fitting software version 2.9 or higher
- CS63 FlexStrip Cable (3-pin)
- Speedlink™, HI-PRO or NOAHlink interface (Speedlink recommended)

ReSound

rediscover hearing

ReSound Live™

Technical Specifications


		LV10-M / LV710-M / LV510-M		LV10-MP / LV710-MP / LV510-MP		
		IEC 118-0 Ear Simulator	IEC 118-7 2cc coupler	IEC 118-0 Ear Simulator	IEC 118-7 2cc coupler	
Reference test gain (60 dB SPL input)	1600 Hz* / HFA**	34	29	41	37	dB
Full-on gain (50 dB SPL input)	Max.	53	43	60	50	dB
	1600 Hz* / HFA**	45	38	53	46	dB
Maximum output (90 dB SPL input)	Max.	119	109	129	118	dB SPL
	1600 Hz* / HFA**	114	106	123	115	dB SPL
Total harmonic distortion	800 Hz	1,2	0,8	1,7	1,0	%
	1600 Hz	1,1	0,9	1,0	0,8	%
Equivalent input noise, w/o Noise reduction		28	26	27	25	dB SPL
1/3 Octave Equivalent Input Noise, w/o Noise reduction		14	-	14	-	dB SPL
Frequency range (DIN 45605*/ANSI**)		130-6800	100-6760	160-6000	100-5990	Hz
Current Drain		0,85	0,88	0,9	0,98	mA
Typical Battery life time (Battery type 10A)		106	102	100	92	hrs

Data in accordance with IEC 60118-0, IEC 60118-7; Supply Voltage 1.3 V.


* IEC 118-0, Ear Simulator

** IEC 118-7, 2cc coupler


Maximum Output (OSPL 90)

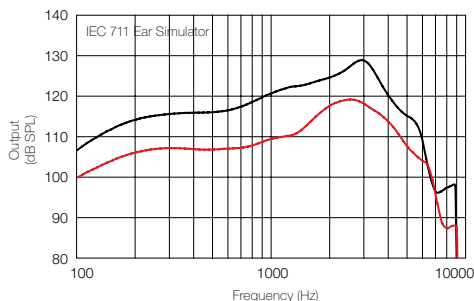


Maximum Output (OSPL 90)

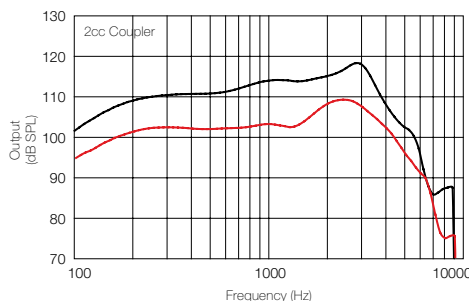


Full-On and Reference Test Gain

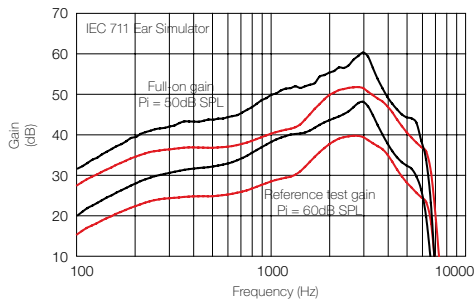

Full-On and Reference Test Gain

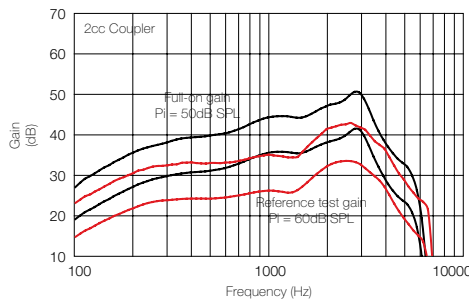

Input/Output Response

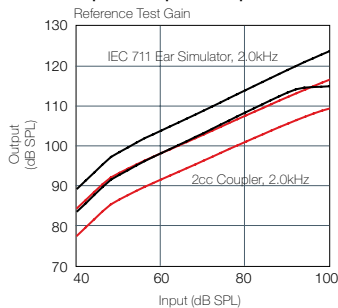

Full-On Gain Parameter Settings*

	250 Hz	500 Hz	750 Hz	1 kHz	1.5 kHz	2 kHz	3 kHz	4 kHz	6 kHz									
G[80]	24	31	24	31	24	31	26	37	26	37	26	35	26	37	24	27	24	27
G[50]	40	46	40	46	40	46	42	52	42	52	42	50	42	52	40	42	40	42

Reference Test Gain Parameter Settings for 118-0*

	250 Hz	500 Hz	750 Hz	1 kHz	1.5 kHz	2 kHz	3 kHz	4 kHz	6 kHz									
G[80]	17	24	17	24	17	24	19	30	19	30	19	28	19	30	17	20	17	20
G[50]	33	39	33	39	33	39	35	45	35	45	35	43	35	45	33	35	33	35

Reference Test Gain Parameter Settings for ANSI and 118-7*

	250 Hz	500 Hz	750 Hz	1 kHz	1.5 kHz	2 kHz	3 kHz	4 kHz	6 kHz									
G[80]	20	28	20	28	20	28	22	34	22	34	22	32	22	34	20	24	20	24
G[50]	36	43	36	43	36	43	38	49	38	49	38	47	38	49	36	39	36	39

*Settings in accordance with Aventa fitting software

ReSound

rediscover hearing

Patents pending.

All specifications are subject to change without notice.

17185800-GB-09.07-Rev.A